




Align Technology iTeroTM



Intraoral Workflow Assistance Guide

800-320-0230 * www.bayoustatelab.com


How to Fill Out an Rx

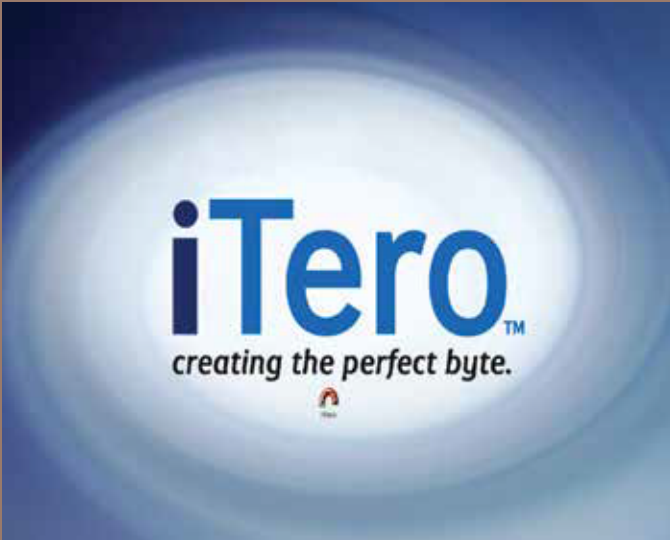

Press the iTero icon



Click "Rx" to create a new case

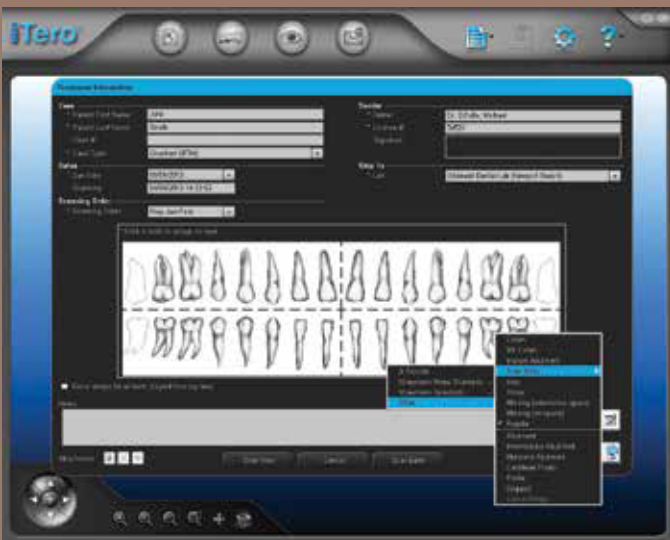

Fill out the required information indicated by *




Fill out notes:

- Implant system
- Implant diameter
- Restoration type


Click "Ok"




How to Prep Scan

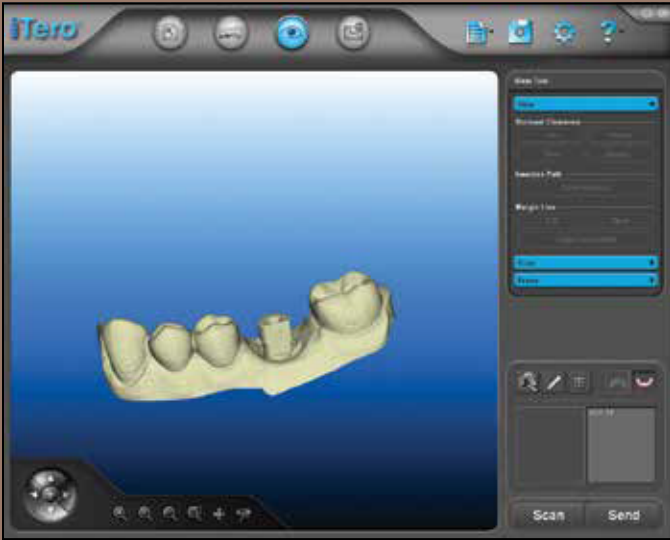


Take prep scan, press "Next"

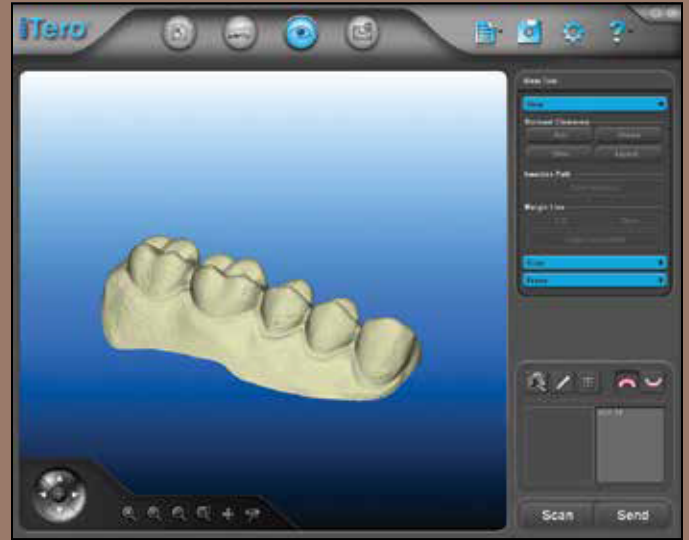

Take opposing arch scan, press "Next"

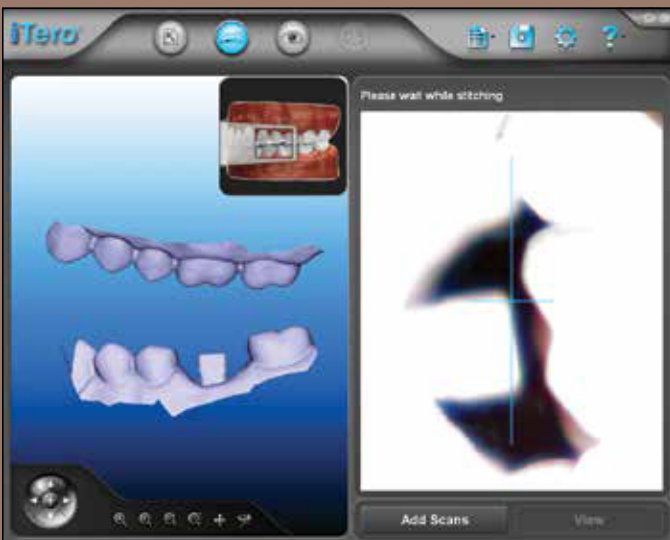


Take bite scan as indicated

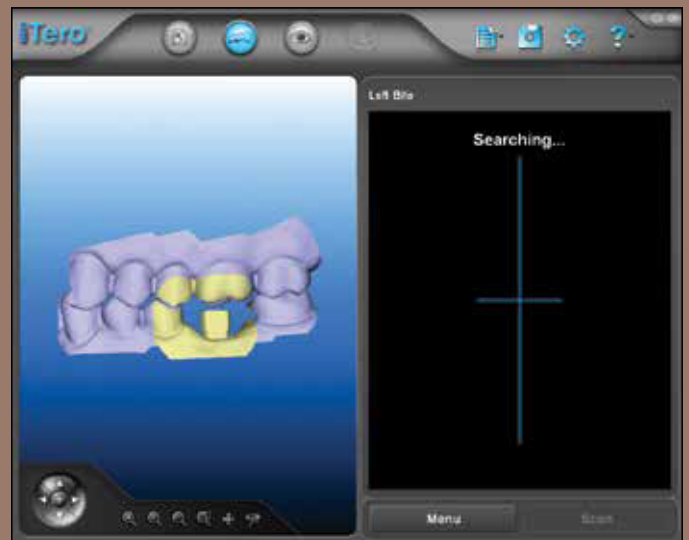


Bite scan complete, press "View"



How to Send Your iTero Files to Bayou State Crown & Bridge Lab Inc.

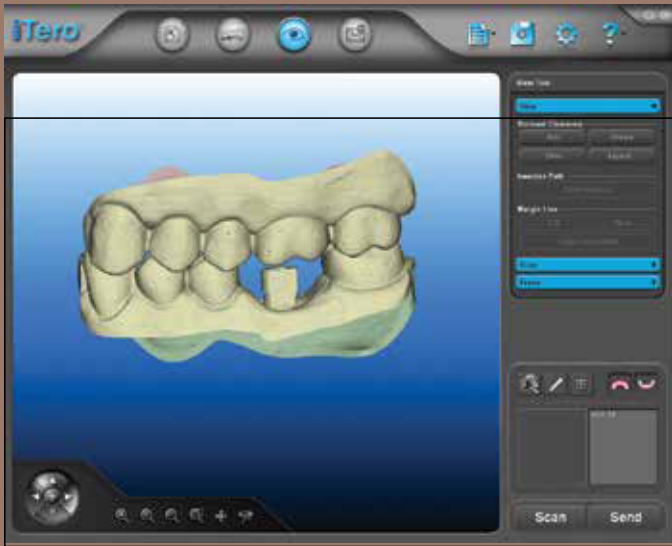

When the digital model is complete, press "Send"

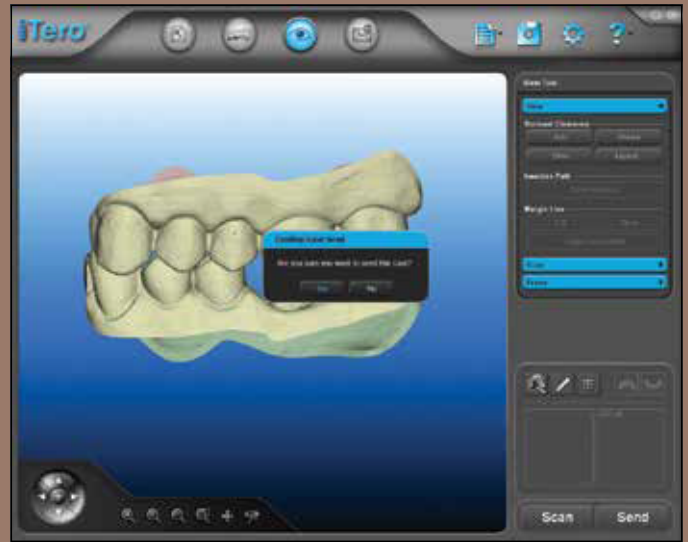

A box will appear to confirm you want to send the case, click "Yes"



800-320-0230 * www.bayoustatelab.com